

Eighteeth. Endo Motor. E-Connect.

E-CONNECT Pro CORDLESS ENDO MOTOR

First Endo Motor with
Five Reciprocation Modes (M0)



Reverse Reciprocation – One Mode
Counter clockwise & clockwise movement
Fwd 30° – Rev 150°
For reciprocating files

Forward Reciprocation – Four Modes
Clockwise & counter clockwise movement

Fwd 150° – Rev 30° | Fwd 180° – Rev 30°
Fwd 210° – Rev 30° | Fwd 250° – Rev 30°
For all continuous rotation files

E-CONNECT Pro CORDLESS ENDO MOTOR

First Endo Motor with
Five Reciprocation Modes (M0)



Reverse Reciprocation – One Mode
Counter clockwise & clockwise movement
Fwd 30° – Rev 150°
For reciprocating files


Forward Reciprocation – Four Modes
Clockwise & counter clockwise movement

Fwd 150° – Rev 30° | Fwd 180° – Rev 30°
Fwd 210° – Rev 30° | Fwd 250° – Rev 30°
For all continuous rotation files

E-CONNECT Pro CORDLESS ENDO MOTOR

Six Editable Memory Programs (M1 to M6)

- 360° 360° Clockwise Rotational Movement
- Editable SPEED / RPM Setting
- Editable Torque Setting
- Forward Direction of Rotation
- Reverse Direction of Rotation

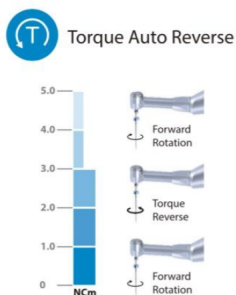

E-CONNECT Pro CORDLESS ENDO MOTOR

Endodontic 1 : 1 Contra Angle Head



- 340° Adjustable Hand Piece
- Detachable Head
- Miniature Head
- 134° C Autoclavable


E-CONNECT Pro CORDLESS ENDO MOTOR

Unique Features of E-Connect Pro Endomotor

- Torque Calibration
- Digital OLED Display
- Low Noise Motor
- Ambidextrous* Display
- Powerful 1600 mAh Lithium-ion Battery



- Integration Option
- Screen Cast / Mirroring*

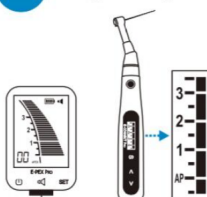

* User can rotate the Endomotor Screen towards left or right

* In connect function Endomotor mirrors Apex Locator's screen

Endodontic Equipment: Eighteeth. Endo Motor. E-Connect.

Eighteeth. Endo Motor. E-Connect.

Rating: Not Rated Yet

[Ask a question about this product](#)